



# USE NQAD IN TANDEM WITH UV TO OBTAIN A COMPLETE PROFILE OF SAMPLE COMPONENTS

Application note # QAN-0077

## PARAMETERS

### NQAD


GAS: Air  
FILTER: 1.3 sec  
GAIN: 1X  
EVAPORATION TEMP: OFF

### HPLC

MOBILE PHASE A: Water  
MOBILE PHASE B: ACN  
METHOD: 0-95% B in 40 minutes

FLOW RATE: 1.0 mL/min  
COLUMN: Alltima amino 100A, 250 x 4.6 mm, 5u  
COLUMN TEMP: Ambient  
INJECTION VOLUME : 10 µL  
ANALYTE: Stash® black chai double spice tea  
CONCENTRATION: 5 g tea steeped for 5 minutes in 100C water (20 mL)

## CHROMATOGRAMS

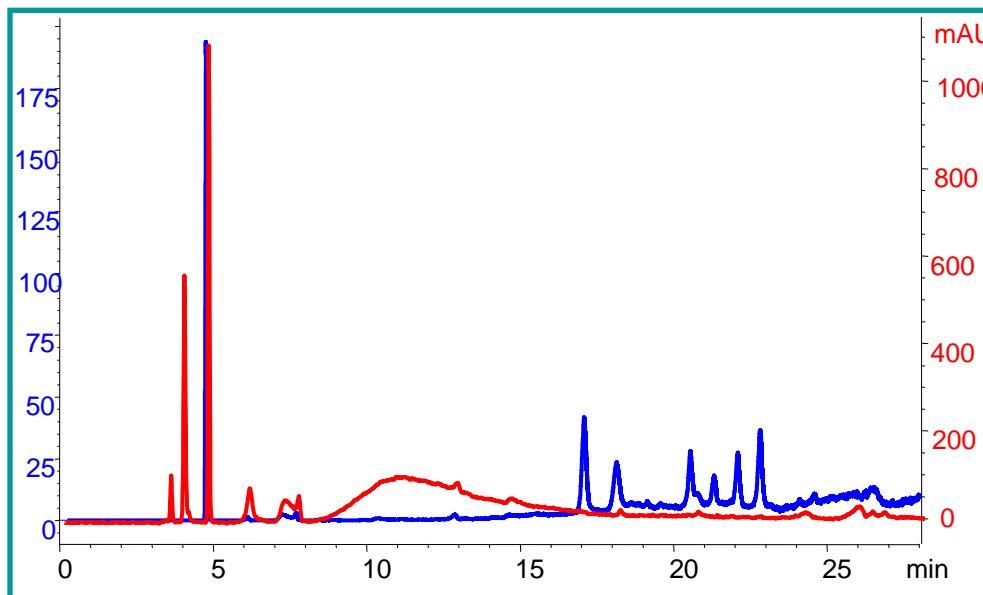

**QUANT TECHNOLOGIES, LLC**  
1463 94<sup>TH</sup> Lane NE  
Blaine, MN, USA 55449  
763-398-0508  
[www.Quant-NQAD.com](http://www.Quant-NQAD.com)

— NQAD  
— UV, 280